**eldom** **linea**

- ELEKTRONICZNA WAGA OSOBOWA •
- ELECTRONIC SCALE •
- ELEKTRONICKÁ OSOBNÍ VÁHA •
- ELEKTRONISCHE PERSONENWAAGE •
- ЭЛЕКТРОННЫЕ ВЕСЫ •
- ELEKTRONICKÁ OSOBNÁ VÁHA •
- ELEKTRONIKUS SZEMÉLYMÉRLEG •
- BÁSCULA ELECTRÓNICA DE BAÑO •

**model: TWO140**

**Eldom Sp. z o.o.**  
Pawła Chromika 5a, 40-238 Katowice, POLAND  
tel: +48 32 2553340 , fax: +48 32 2530412

**eldom.eu**

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## SAFETY

Please read all these instructions carefully before using this appliance.

1. This device may be used by children of 8 years of age or more and by persons with limited physical and mental capacity, as well as by persons without experience and knowledge of the device, provided that supervision or instruction is ensured as regards operation of the equipment in a safe way, so as to make the related hazards understandable. Children should not play with the device. Children without supervision should not clean and maintain the appliance.
2. The device is intended only for household use.
3. Do not stand on the scales when it is wet.
4. Do not stand on the scales with wet feet.
5. Store the scales in a place shady and dry place.
6. When in storage, the scales should not be placed under other items.
7. The accuracy of reading of the scales depends on the way it is used. Use the scales carefully. Do not jump on the scales, do not drop the scales.
8. Use damp, soft cloth without chemical agents to clean, then wipe dry. Do not immerse in water.
9. Only authorised service centres may repair the device. Any modernisations or using of non-original spare parts or elements of the device is forbidden and threatens safety of its use.
10. Eldom Sp. z o.o. shall not be liable for any damages resulting from improper use of the device.

### TECHNICAL DATA

- weight capacity 150kg (330lb)
- resolution 100g (0,2lb)
- accuracy of measurement: 1kg
- power: CR2032

This product is specially designed to test body fat/hydration, help keep a close eye on the daily change of your body fat, hydration and weight, providing information for your reference. Base on the analysis of your height, weight and age etc., it tells your body fat status and gives you an idea of your body condition. Please read this instruction carefully before use.

### During Measuring


1. Always remove your shoes and socks, wipe clean your feet before using.
2. It is advisable to take measurement at the same time of the day.
3. Measured result could be misleading after intensive exercise, excessive dieting or under extreme dehydration condition.
4. Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
5. The body fat/hydration measured data of the following people may have deviation:

- Children under 10 or adult over 70 (can only use the device under normal weighing mode)
  - Body building athletes or other occupational athletes.
6. People with the following conditions should not use the device:
    - Pregnant women
    - Those on dialysis treatment
    - **Those who use heart pacemaker or with other implanted medical device.**

### Safety Warning

The Fat% Hydration% Muscle% content evaluated by this device is only for reference (not for medical purpose). If your Fat% Hydration% content is over/ below normal level, please consult your doctor for more advice

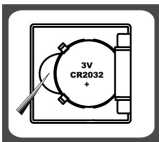
### BATTERY


Used batteries should be discarded in special collection bins. Please remove the isolation sheet from the battery contact before use.

### BATTERY REPLACEMENT:

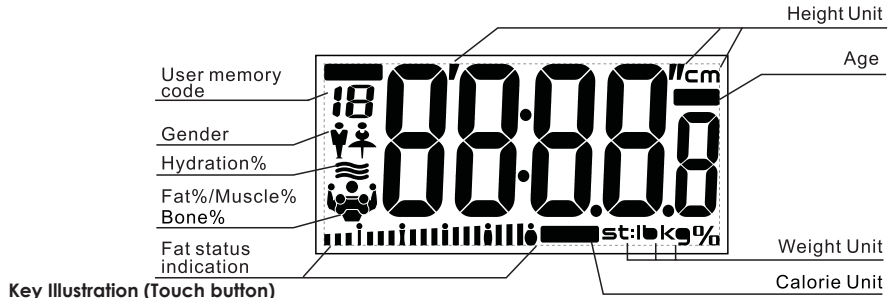1. Open the battery cover at the bottom of the device.
2. Carefully lever the battery and take it out.
3. Insert a new battery and press it gently until the retaining clip snaps in place. Please pay attention to the polarization.
4. Close the lid.

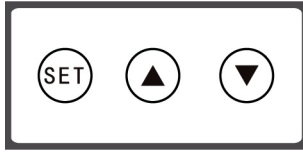Key Illustration (Touch button)

-  Turn on the scale and prompt to next parameter. (Please refer to the Setting Personal Parameter overleaf)
-  Advance the value by one interval each press. Continuous pressing will speed up the advance of value.
-  Decrease the value by one interval each press. Continuous pressing will speed up the decrease of value.

**TO SET WEIGHT UNIT**

- Press "SET" and it will get ready for use when screen shows "0,0kg" after "8888" and you can press ▲ and ▼ to switch within different units (kg/lb)

**FOR WEIGHT ONLY**

- Press "SET"
- Step on the scale gently, position your feet evenly near the center and keep your body still until the scale stabilizes your weight.
- The scale will be auto off if no further using.

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-  Press ▲ or ▼ key to select a number from P0 to P9 as your personal code. Please press "SET" to confirm
-  Press ▲ or ▼ key to choose  Male or  Female. Please press "SET" to confirm
-  Press ▲ or ▼ key to reach your target age (10-80). Please press "SET" to confirm.
-  Press ▲ or ▼ key to reach your target height (100-240 cm). Please press "SET" to confirm.
-  Please press "SET" to confirm.

**Body fat and hydration ranges table.**

Age	Female		Male		Display
	Fat %	Hydration %	Fat %	Hydration %	
<=30	4.0-20.5	66.0-54.7	4.0-15.5	66.0-58.1	Underfat
	20.6-25.0	54.6-51.6	15.6-20.0	58.0-55.0	Healthy
	25.1-30.5	51.5-47.8	20.1-24.5	54.9-51.9	Overfat
	30.6-80.0	47.7-13.8	24.6-80.0	51.8-13.8	Obese
>30	4.0-25.0	66.0-51.6	4.0-19.5	66.0-55.4	Underfat
	25.1-30.0	51.5-48.1	19.6-24.0	55.3-52.3	Healthy
	30.1-35.0	48.0-44.7	24.1-28.5	52.2-49.2	Overfat
	35.1-80.0	44.6-13.8	28.6-80.0	49.1-13.8	Obese

Step on the scale gently and keep your bare feet in touch with the electrodes. Stand still and erect to get your weight measured and your body fat, hydration, muscle and bone percentage will be displayed respectively on the screen within seconds. These results will be repeated for two times and the scale will be auto off no further using.

The personal data will be stored once you set. Just confirm your code when you use scale again and your other information will be showed so there is no need to program them again.

### **CLEANING AND MAINTENANCE**

- The appliance cannot be submersed in water
- Do not use any cleansing agents which are strong and damage the surface
- Wipe the housing with a moistened cloth

### **PROTECTION**

- the appliance is made of materials which can be recovered or recycled
- the appliance should be handed over to the appropriate centre of collection and recycling of electrical and electronic equipment

### Body muscle and bone ranges table

Age	Female		Male	
	Muscle %	Bone %	Muscle %	Bone %
16-50	>31	>8	>38	>10

### BMI parameter

Low	Normal	Overweight	Fat
<18.5	18.5~25	25~30	>30

### Warning Indications

Low Battery Indication.  
The battery power is running low, please replace with a new battery.

Over-load Indication  
The weighing subject on the platform exceeds the maximum capability of scale. Please step off to avoid damage.

Rate too low level of fat.

Too high a level of fat.

Message redo the measurements.

If body percentage is less than 5%, then screen will show "FATL". If more than 50% it will display "FATH" and "ERR2" indicates that the measuring makes an error.

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