
Oral-B®

KIDS

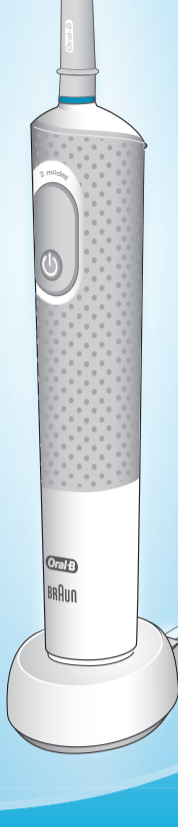powered by **BRAUN**

Charger
Type 3757
Handle
Type 3710

90714170

- | | |
|----|----------------------------|
| UK | 0 800 731 1792 |
| IE | 1 800 509 448 |
| PL | 801 127 286
801 1 BRAUN |
| CZ | 800 11 33 22 |
| SK | 800 333 233 |
| HU | (06-1) 451-1256 |
| HR | 01 / 66 90 330 |
| SI | 080 2822 |
| LT | (8 5) 205 1272 |
| LV | 67798667 |
| EE | 667 5047 |
| IL | 1-800-666-775 |

Internet:










www.oralb.com

www.braun.com

www.service.braun.com

Braun GmbH
Frankfurter Straße 145
61476 Kronberg/Germany/
Germania/ Германия/ Немацка

EN/PL/CZ/SK/HU/HR/SI/LT/LV/EE/IL
90714170/X-19

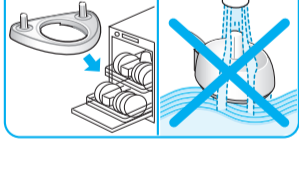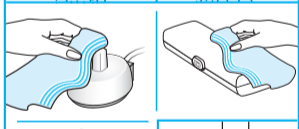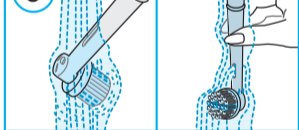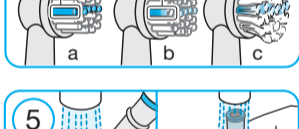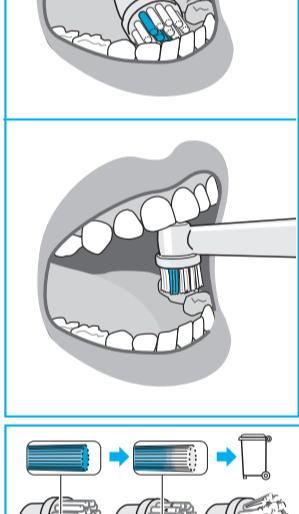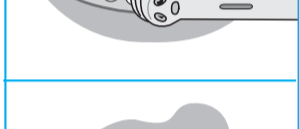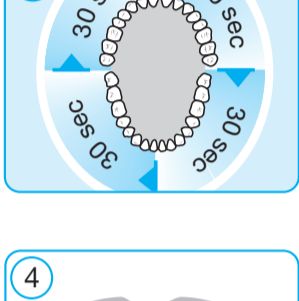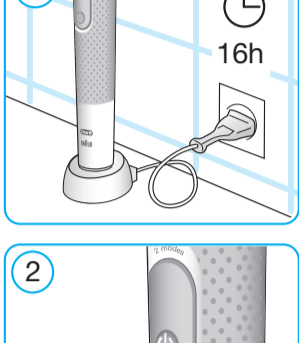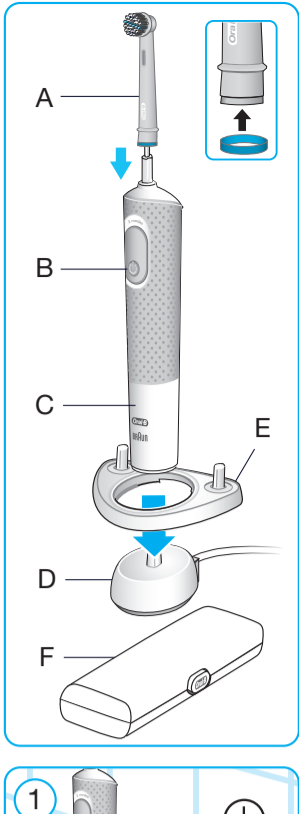English

Dear parents,
Please read the usage instructions carefully before coaching your child how to use this toothbrush. We recommend that parents supervise their children when they start using the appliance to help them to become familiar with it.

IMPORTANT

- Periodically check the entire product/cord for damage. A damaged or non-functioning unit should no longer be used. If the product/cord is damaged, take it to an Oral-B Service Centre. Do not modify or repair the product. This may cause fire, electric shock or injury. Never insert any object into any opening of the appliance.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the

appliance in a safe way and understand the hazards involved.

- For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer. Only use charger provided with the appliance.

WARNING

- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- This appliance contains batteries that are non-replaceable. Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Follow correct instructions when applying stickers on handle. Ensure valve on back of the handle is not covered.

To avoid brush head breakage which may generate a small parts choking hazard or damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fit properly. Never use without a brush head.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible
- Replace brush head every 3 months or sooner if brush head becomes worn.
- Clean brush head properly after each use (see «Cleaning recommendations» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.

Description

- A Brush head
- B On/off switch (mode selection)
- C Handle
- D Charging unit

Accessories (depending on model):

- E Brush head holder
- F Travel case (design may vary)

Note: Content may vary based on the model purchased.

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Noise level: ≤ 65 dB (A)

Charging and Operating

This toothbrush is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet and place the toothbrush handle (C) on the charging unit. A full charge takes about 16 hours and enables up to eight days of regular brushing (twice a day, 2 minutes) (1).
- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Using the toothbrush

- Wet brush head and apply kids' toothpaste. Guide the brush head to the teeth; then turn on the handle by pressing the on/off switch (B) (2).
- Start with the lower teeth: Slowly guide the brush head from tooth to tooth; first the outsides, then the insides, finally the chewing surfaces (4). Do not press too hard, simply let the brush do all the work.
- A short stuttering sound every 30 seconds reminds to brush all four quadrants of the mouth equally. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time (3).
- If the toothbrush is turned off during brushing, the elapsed brushing time will be memorised for 30 seconds. When pausing longer than 30 seconds the timer resets.
- Turn off the handle while the brush head is still in the mouth.

In the first few days of using any electric toothbrush, gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or your dental hygienist.

Brushing modes

«**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas

«**Daily Clean**» – Standard for daily mouth cleaning

How to switch between modes:

The toothbrush automatically starts in the «Sensitive» mode.

To switch to other mode, press the on/off button (B). Switch off the toothbrush by holding the on/off button until the motor stops.

Replacement brush heads: Oral-B Stages Power

Available at your retailer or at an Oral-B Braun Service Centre.

The brush head is provided with a row of blue INDICATOR® bristles (4a) helping you to monitor brush head replacement need. With proper brushing twice a day for two minutes, the colour will fade halfway in a 3-month period (4b). If bristles begin to splay before the colour starts to recede, this indicates that your child exerts too much pressure on teeth and gums (4c).

Cleaning recommendations

After use, rinse the brush head for several seconds under running water with the handle switched on. Switch

off and remove the brush head. Rinse both parts separately running water and wipe them dry. Charging unit (D) and travel case (F) should be cleaned with a damp cloth only. Never place charging unit in water. Store the clean and dry toothbrush/brush heads in the travel case. Brush head holder (E) is dishwasher-safe (5).

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. Depending on availability replacement of a unit might result in a different color or an equivalent model. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Oral-B Service Centre.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

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