# SHONA SONIC PRO

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
**Výkon:** 1.8W

## SHONA SONIC PRO TOOTHBRUSH

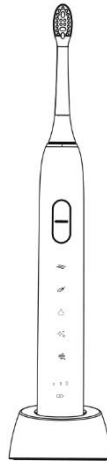Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and maintenance shall not be made by children without supervision.

Children shall not play with the appliance.

The supply cord cannot be replaced. If the cord is damaged, the appliance should be scrapped.

The appliance is only to be used with the power supply unit provided with the appliance.

### **S710 SONIC ELECTRIC TOOTHBURSH IMPORTANT SAFEGUARDS**

Shona Sonic Pro has been carefully designed to offer you and your family a new and unique brushing experience, which is both safe and effective. However, when using the electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

**Read all instructions before using the appliance**

Reduce the risk of electrocution:

1. Always unplug this product immediately after using.
2. Do not use while bathing
3. Do not place or store the product where it can fall or pulled onto a tub or sink
4. Do not place in or drop into water or other liquid
5. Do not reach for the product that has fallen into water unplugged immediately

### **WARNING**

To reduce the risk of fire, electric shock, or injury:

Do not leave appliance when plugged in. Unplug from outlet when not in use and before servicing

Do not allow to be used as a toy. Close attention is necessary when used by or near children

Use only as described in this manual. Use only manufacturer's recommended attachments

Do not use with damaged cord or plug. If appliance is not working as it has been dropped, damaged, left outdoors, or dropped into water, returning it to a service center

Do not pull or carry by cord, use cord as a handle, close a door on cord, or pull cord around sharp edges or corners

Do not run appliance over cord. Keep cord away from heated surfaces

Do not unplug by pulling on cord. To unplug, grasp the plug, not the cord

Do not handle plug or appliance with wet hands

Turn off all controls before unplugging

Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where air openings may be blocked. Keep the air openings free of lint, hair, and the like.

Never use while sleeping or drowsy

Never drop or insert any object into any opening or hose

Do not use outdoors or operate where aerosol products are being used or where oxygen is being administered

Connect this product to a properly grounded outlet only. See grounding instructions

Unplug this product before filling. Do not overfill.

**SAVE THESE INSTRUCTIONS HOUSEHOLD USE ONLY PLEASE READ BEFORE USE:**

This product has no user serviceable parts

Store the charging unit onto a dry and safe place when it is not in use


Untwisted the power cord before use

In case of the product failing to function, you should not dismantle, exchange or manipulate any part of it, otherwise there is risk of an electric shock. You should take it to the seller.

If the supply cord is damaged, it must be replaced by the manufacturer in order to avoid a hazard. You should contact the seller

If the product is dropped, the brush head should be replaced before the next use even if no damaged is visible

**LIST OF COMPONENTS**

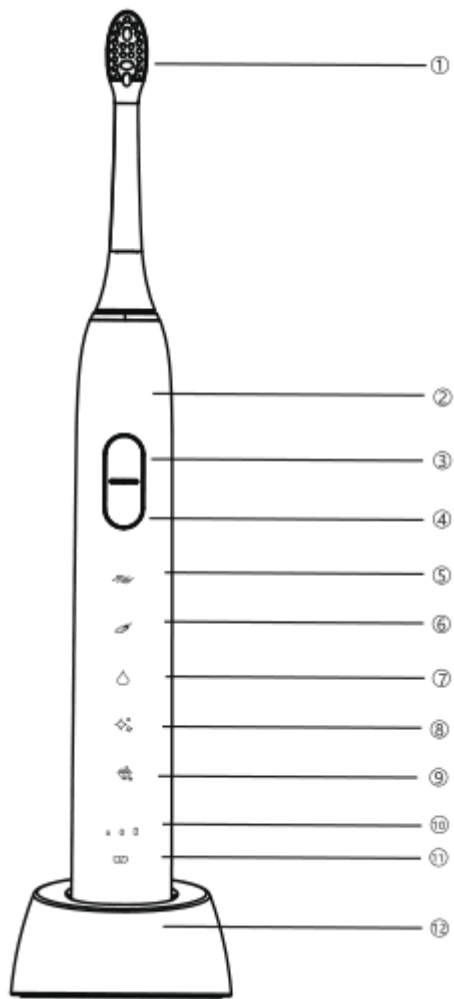1. Brush Head
2. Toothbrush Handle
3. Power on/off
4. Brushing Strength Button
5. Clean
6. Sensitive
7. Massage
8. White
9. Polish

10. Brushing Strength Indicator

11. Charging Indicator

12. Charging Base

## **CONNECTING AND CHARGING**

Shona Sonic Pro has a waterproof handle. It is electrically safe and is designed for use in the bathroom. But it's not meant to be used while underneath the water.

Plug the adaptor into an electrical outlet and put the brush handle onto the charging base, the charging indicator should be lit.

The first full charge needs 16 hours and later full charge needs 12 hours. That enables 30 days brushing. Power indicator light turns red and flashing means battery low power. indicator light keeps flashing red when charging, and fully charged when indicator turns green.

Pull plug out of outlet for saving energy, or leave it and keep the brush on charging base, there is no risk of overcharging of battery.

If not use for a long time and when you charge it, the power indicator will not be lit immediately which is normal. After pre-charging 30 minutes, the power indicator will be lit.

## **USING THE Shona Sonic Pro**

Attached a brush head to the handle. You can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on.

## **BRUSHING TECHNIQUE:**

Guide the brush head slowly from tooth to tooth

Hold the brush head in place for a few seconds before moving on to the next tooth.

Brush the gums as well as the teeth. Start with outside of the teeth, then the insides ending with the chewing surfaces.

Do not scrub or press hard, simply let the brush do all work.

## **BRUSHING MODES**

Shona Sonic Pro offers different brushing modes and strength for varying oral care needs:

Clean with 3 different brushing strengths Brush head in high frequency, remove bacteria on teeth.

Sensitive with 3 different brushing strengths- For sensitive teeth and gums, brushing your teeth in gentle frequency.

Massage with 3 different brushing strengths-Brush teeth in extremely gentle frequency.

White with 3 different brushing strengths- Clean and Polish mode switching continuously and whitening teeth.

Polish with 3 different brushing strengths- High strength and high frequency, polishing and brightening teeth.

## **TWO MINUTES AUTO-TIMER:**

Product is designed to reach the proven and best brushing timing. After two minutes of brushing, the device stops automatically.

## **HOW TO SWITCH BETWEEN BRUSHING MODES:**

Power on/off and brushing mode button. Keep pressing this button within 5s, the operating sequence is power on, in clean mode, sensitive mode, massage mode, white mode, polish mode, end power off. Keep pressing the button 5s over, it will power off.

The S710 automatically starts in the clean mode when pressing the brushing mode button once

To shift from one mode to another mode, press the brushing mode button a second time

When you press the brushing mode button, the brushing mode indicator will be light on, indicating the selected working mode

2 minutes in Clean, Sensitive, Massage, White and Polish modes and 30s interval timer guide.

#### **CLEANING RECOMMENDATIONS:**

After brushing, rinse your brush head thoroughly under running water

Remove the brush head from the handle and clean both parts separately under running water.

Then wipe them dry before you put the handle onto the charging unit

Before cleaning your charging unit, make sure you unplug the adaptor. Only use a dry or damp-dry cloth to wipe the surface of the charging unit

Place the charging unit onto a dry and safe place

Do not plug in the charging unit until it is completely dry

#### **DISPOSAL**

Important: do not clean the charging unit under running water or immerse it into water. Place the charge station on dry and safe place

This product contains rechargeable battery. In the interest of protection the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposable can take place at appropriate collection points provided in your country.

## PARAMETER

Rated input: DC5V 1A

Main working voltage: 3.0V-4.2V

Battery capacity:800mAh

Rated Power: 1.8W