# BEAUTIFLY

B-Scrub Slim

# PRO

series

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
Multifunctional skin scrubber

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User manual

## Opis części

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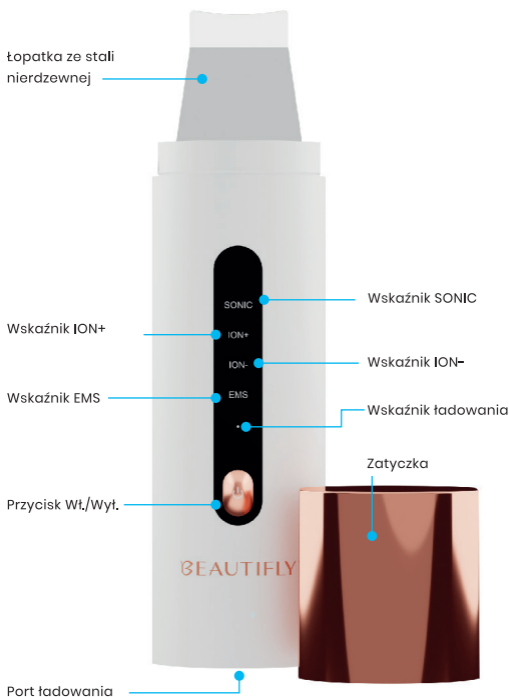
## Part description
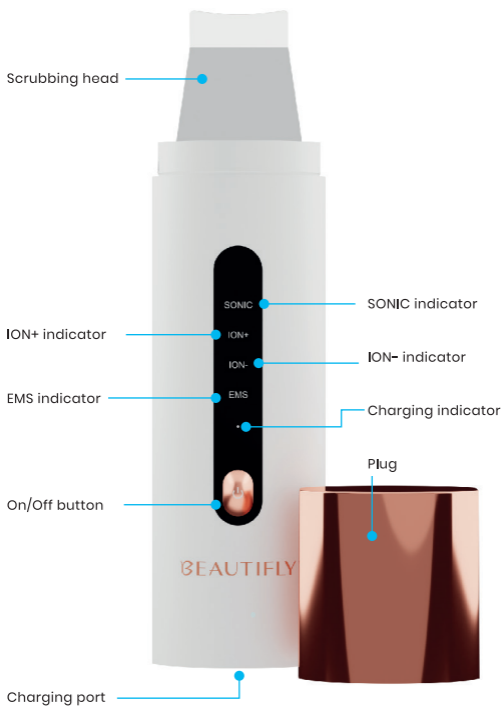
### Technical data B—Scrub PRO Slim

<b>product dimension:</b>	170*44*20mm
<b>net weight:</b>	~122g
<b>lithium battery:</b>	3.7V === 700mAh
<b>input power:</b>	5W
<b>charging voltage:</b>	5V === 1000mA
<b>working temperature:</b>	+5°C~+40°C
<b>relative humidity:</b>	≤80%

## Function

Multifunctional cavitation peeling B—Scrub PRO Slim for skin care. It combines the technology of ultrasound and positive ions for atomization and cleaning of dead, calloused epidermis. The combination of ultrasounds and negative ions promotes the absorption of the essence contained in used creams, ampoules or serums. Ultrasound and EMS stimulate the production of collagen and elastin, responsible for skin tightness. Use 1-2 times a week to thoroughly cleanse the skin and improve the effectiveness of daily care.

### Ultrasonic cleansing (SONIC)

When the skin is moistened with water, high-frequency ultrasonic vibrations atomize the water, completely softening the dead skin.

### Cleansing (SONIC & ION +)

After ultrasonic cleansing in the SONIC mode, the skin should still be moist. Positive ions help to export contaminants deep in the pores of the skin.

### Moisturizing (SONIC & ION-)

Attracting anions use the microelasticity of human skin, using a peeling spatula, they effectively import nutrients to the surface of the base layer of the skin, leaving it nourished and moisturized.

### Rejuvenation (SONIC + EMS)

Microcurrents generated by ultrasound and EMS stimulate cell activity. This process stimulates the skin to produce collagen and elastin. As a result, the skin becomes firm and smooth.

## Description of modes

No.	Status	Indicator
1.	Cleansing	Turn on the SONIC device, the peeling blade vibrates at high frequency
2.	Cleansing	Turn on the device and press the on / off / mode button on ION + 2 times, the peeling paddle vibrates at high frequency
3.	Moisturizing	Turn on the device and press the on / off / mode button on ION- 3 times, the peeling paddle vibrates at high frequency
4.	Anti-age	Turn on the device and press the on / off / mode button 4 times in EMS mode, the peeling head vibrates in high frequency
5.	Auto-off	Beeps 2 times, lights off, device off.
6.	Manual off	Press on/off/mode to skip all modes till turn it off.
7.	Low battery	Beeps 3 times, button turns red and flash.
8.	In charging	On/off/mode button turns red and keep flashing.
9.	Full charged	On/off/mode button turns green.

## Usage method

### Cleansing the epidermis (Ultrasounds)

1. Before starting the treatment, thoroughly cleanse the skin and remove makeup.
2. Then moisten the skin of the face with a cotton swab or an atomizer. For this purpose, you can use saline, tonic, water or hydrolate. During the treatment, the skin must be moist all the time. The treatment cannot be performed on dry skin.
3. Turn on the device and start the SONIC mode.
4. Gently move the spatula from bottom to top, lengthwise facial contours at a 45-degree angle. Repeat the action carefully, several times over the entire face.

### Cleansing the epidermis with ionization (Ultrasound + Ionization +)

1. Before starting the treatment, remember to keep the skin moist all the time. For this purpose, you can use saline, tonic, water or hydrolate. The treatment cannot be performed on dry skin.
2. Turn on the device and enter ION + mode.
3. Gently move the spatula from bottom to top, lengthwise facial contours at a 45-degree angle. Repeat the action carefully, several times over the entire face.
4. Wipe your face with a tonic.

### The introduction of nutrients with ultrasound (Ultrasound + Ionization-)

1. After cleansing the skin of the face, apply a skin care product (cream, serum, ampoule, gel).
2. Turn on the device and start the ION mode.
3. Gently move the spatula from the bottom up along the contours of the face. The spatula should be placed in the opposite direction to that in step one. Repeat the action carefully, several times over the entire face.

If you have any doubts how to proper hold the device, please visit our chan on the YouTube platform.

 /beautifyeu

### EMS (Ultrasound + EMS currents)

1. After caring for your facial skin, turn on the device and enter the EMS mode.
2. To improve skin elasticity, move the spatula in an upward motion along the contours of your face at a 45-degree angle. Repeat the action carefully, several times over the entire face.
3. The device will turn off automatically after 10 minutes of operation
4. After the device is finished, wipe it with a cotton cloth damp cloth.

## Cautions

1. Person with electronic human organ are forbidden to use it.
2. Person with skin surgery within 1 month or skin sunburn is forbidden to use.
3. Forbidden to use by person with inflammation of the skin, scar, damaged and electromagnetic sensitive is forbidden to use.
4. Take off watch, glasses and other metal objects before operation.
5. Avoid shovel head close to eyeball.
6. Don't operate in the same area for long time.
7. Do not use the device for pregnant women.

## Cleaning

Rinse the shovel head and wipe it with a dry soft cloth after using. Do not immerse the device into water for long time.

## Storage

1. Keep it dry and clean.
2. Keep it away from high temperature.
3. Keep it out of reach of children.

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[/beautifyyou](https://www.youtube.com/beautifyyou)

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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### English

Hereby, Beautifly Sp. z o. o. declares that the device equipment Beautifly B-Scrub PRO Slim is in compliance with Directive 2014/30/UE.

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## TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
3. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
4. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address [serwis@rgmedia.pl](mailto:serwis@rgmedia.pl)

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