DE Anti-Aging Gesichtspflege Gebrauchsanweisung	2	IT Cura del viso anti-aging Istruzioni per l'uso.....	34
EN Pureo Ionic Skin Care Instructions for use	10	TR Yaşlanma önleyici yüz bakımı Kullanım kılavuzu	42
FR Soin du visage anti-âge Mode d'emploi.....	18	RU Прибор для ухода за кожей лица Инструкция по применению.....	50
ES Cuidado facial antiedad Instrucciones para el uso.....	26	PL Pielęgnacja twarzy Anti-Aging Instrukcja obsługi.....	57

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Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Warning

- The device is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of eight and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- The device must not be used by persons with medical implants (e.g. heart pacemakers) Otherwise their function could be impaired.
- Do not use the device if it shows signs of damage or does not function properly. In these cases, contact Customer Services.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.

Dear customer,

Thank you for choosing a product from our range. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage, beauty, air and baby.

With kind regards
Your Beurer team

Getting to know your instrument

Thanks to its multifunctional ionic technology, the FC 90 Pureo Ionic Skin Care is ideal for cleansing and anti-aging facial treatments. The treatment head can be charged with either positive or negative ions. Electrochemical reactions cause these charged ions to get under the skin's surface into the deeper layers of the skin or to move from the deeper layers of the skin to the skin's surface.

Desincrustation cleansing:

The negatively charged components in the cleansing milk are channelled deep into the skin with the aid of the negative ions. There, they loosen impurities, which are then drawn to the skin's surface by the positive charge of the treatment head.

Anti-aging by iontophoresis:

If the treatment head is negatively charged, the negatively charged anti-aging agents in the Advanced Hydro Lift Cream are pushed away from the treatment head and carried deep into the skin. There, the anti-aging agents have the best possible effect.

Cool and heat function:

The FC 90 Pureo Ionic Skin Care also has a cool and heat function. The heat function relaxes the skin on the face and promotes circulation. It leaves your skin feeling smooth. With the cool function, swellings on the face can be reduced, pores can be tightened and the skin elasticity can be improved.

Contents

1. Included in delivery.....	11	8. Cleaning and maintenance.....	15
2. Signs and symbols.....	11	9. Accessories and replacement parts.....	15
3. Intended use.....	12	10. What if there are problems?.....	15
4. Warnings and safety notes.....	12	11. Disposal.....	15
5. Device description.....	13	12. Technical specifications.....	16
6. Initial use.....	13	13. Warranty / Service.....	16
7. Usage.....	13		

1. Included in delivery

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 FC 90 Pureo Ionic Skin Care
- 1 Mains adapter
- 1 Intensive Cleansing Milk (30 ml) *(You can reorder the cleansing milk specially designed for use with the FC 90 from our Customer Services by quoting the order number 163.394.)*
- 1 Advanced Hydro Lift Cream (15 ml) *(You can reorder the Advanced Hydro Lift Cream specially designed for use with the FC 90 from our Customer Services by quoting the order number 163.392.)*
- 1 Instructions for use

2. Signs and symbols

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	Warning Warning instruction indicating a risk of injury or damage to health		Read the instructions
	Product information Note on important information		Dispose of packaging in an environmentally friendly manner
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE		CE labelling This product satisfies the requirements of the applicable European and national directives.
	Only use indoors		The products demonstrably meet the requirements of the Technical Regulations of the EAEU.
	Manufacturer		Energy efficiency class 6
	Polarity		

3. Intended use

The FC 90 Pureo Ionic Skin Care must only be used to treat human facial skin.

The device is only intended for the purpose described in these instructions for use. Use the device only according to the methods specified. The manufacturer is not liable for damage resulting from improper or careless use.

This device is intended for domestic use only, not for commercial purposes.

4. Warnings and safety notes

Warning

Before using the FC 90 Pureo Ionic Skin Care, read the safety notes through carefully. Non-observance of the following information may result in personal injury or material damage. Store these instructions for use and make them accessible to other users. Make sure you include these instructions for use when handing over the device to third parties.

- If, after using the device, you experience a worsening of the skin, undesired side effects or prolonged redness, stop using the device and consult your doctor.
- Do NOT use the device on acutely or chronically diseased (injured or irritated) skin (e.g. inflamed skin – whether painful or not, reddened skin, rashes, e.g. allergies, burns, bruises, swellings, as well as open and healing wounds).
- Keep the device away from water.
- Only use with the supplied mains adapter.

Electric shock

Warning!

- As with every electrical device, the FC 90 Pureo Ionic Skin Care must be used carefully and cautiously in order to avoid hazards due to electric shock.
- Only use the device with the mains adapter supplied and at the mains voltage specified on the type plate.
- Only plug the mains adapter into the socket when the device is switched off.
- Never use the device if it or any of its accessories are visibly damaged.
- Do not use the device if it has sustained a fall or if water has penetrated the device interior.
- Do not use the device during a storm.
- In the event of defects or operational faults, switch off the device immediately. Pull the plug out of the mains socket.
- Do not pull on the mains cable or on the device to disconnect the mains adapter from the socket.
- Never hold or carry the device by the mains cable.
- Do not clamp or pinch the mains cable.
- Keep the mains adapter away from hot surfaces.
- Ensure that the mains cable and mains adapter do not come into contact with water, steam or other liquids.
- Never reach for a device that has been submerged in water. Pull the mains adapter out of the mains socket immediately.

Repairs

Warning!

- Do not open the device.
- Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. For repairs, please contact Customer Services or an authorised retailer.

Risk of fire

Warning

There is a risk of fire if the device is used improperly or these instructions for use are ignored.

- Never use the device underneath a cover, such as a blanket or pillow.
- Never use the device near highly flammable substances.
- Make sure that the vents are not covered.
- Never store the device while it is still switched on. Always switch the device off before storing it.
- Do not expose the device to high temperatures.

5. Device description

1. Treatment head
2. LCD display
3. ON/OFF button
4. Negative ions button
5. Positive ions button
6. Heat button
7. Cool button

6. Initial use

Warning

Keep packaging material away from children! There is a choking hazard.

- Remove the device from the packaging and check that the device and all accessories are included and not damaged.
- Plug the mains cable into the bottom of the device.
- Ensure that the mains cable does not pose a trip hazard.
- Plug the mains adapter into the socket.

7. Usage

Warning

Before starting the treatment, carry out a skin sensitivity test to check whether your skin is sensitive to the cleansing milk and the Advanced Hydro Lift Cream. To do this, proceed as follows:

- Clean a patch of skin in the crook of your arm thoroughly, apply a little of the cleansing milk and Advanced Hydro Lift Cream to this cleaned patch of skin and allow the product to soak in and dry.
- If the skin shows no reactions within the next 48 hours, you can use the tested products without hesitation.
- If the skin appears to react before the end of the 48-hour period, you should not use the tested products for facial treatment.
- Do not carry out the skin sensitivity test around the eye area.

You can use the device to care for your skin as well as to cleanse it.

During the treatment, hold the device so that your hand covers both chrome-coloured stripes on the side of the device. If you do not do this, the electric circuit will be open and the ionic function will not be possible.

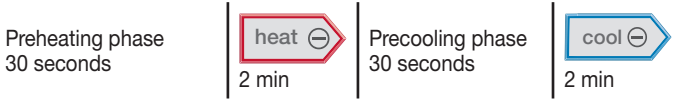
Treatment applications

To ensure the best possible effect, we recommend that you only use the products provided when using the cleansing treatment function and anti-aging treatment function, as these have been perfectly adapted for use with the FC 90.

We recommend using the care treatment function in the morning and a combination of the cleansing and care treatment functions in the evening.


7.1 Care treatment function (mornings):

- Apply the Advanced Hydro Lift Cream to the area of your face that you wish to treat.
- Press the ON/OFF button on the device.
- Press the heat button first and then press the negative ions button . The device is now in the pre-heating phase for 30 seconds. "WAITING" is shown in the display. As soon as the "WAITING" display turns off, the device is ready to use.
- Now use the treatment head to gently massage your face in circular motions.
- After two minutes of treatment, press the cool button and the negative ions button . The device is now in the precooling phase for 30 seconds. "WAITING" is shown in the display. As soon as the "WAITING" display turns off, the device is ready to use. Continue with the treatment for another two minutes.
- To switch the device off again, press the ON/OFF button .

7.2 Combination of cleansing and care treatment functions (evenings):

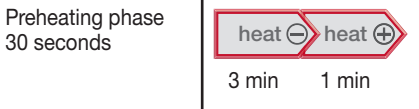
Cleansing:

- Apply the cleansing milk to your face.
- Press the ON/OFF button on the device.
- Press the heat button first and then press the negative ions button . The device is now in the pre-heating phase for 30 seconds. "WAITING" is shown in the display. As soon as the "WAITING" display turns off, the device is ready to use.
- Now use the treatment head to gently massage your face in circular motions.
- After three minutes, press the positive ions button and continue the treatment for another minute. After the cleansing treatment, wash your face with lukewarm water.

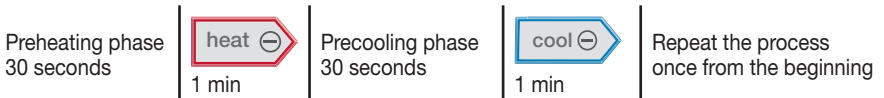Care:

After you have completed the cleansing treatment function, continue with the care treatment as follows:

- Clean the treatment head with a lint-free cloth until there is no trace of cleansing milk left on the head.
- Apply the Advanced Hydro Lift Cream to the area of your face that you wish to treat.
- Press the heat button first and then press the negative ions button . The device is now in the pre-heating phase for 30 seconds. "WAITING" is shown in the display. As soon as the "WAITING" display turns off, the device is ready to use.
- Now use the treatment head to gently massage your face in circular motions.

- After one minute of treatment, press the cool button . The device is now in the precooling phase for 30 seconds. "WAITING" is shown in the display. As soon as the "WAITING" display turns off, the device is ready to use. Continue with the treatment for another minute.
- Repeat this heat and cool treatment as described above for one minute each. If necessary, apply a little more Advanced Hydro Lift Cream.
- To switch the device off again, press the ON/OFF button .

8. Cleaning and maintenance

- When you are not using the device, store it in a cool, dry place away from direct light.
- Carefully clean the device after every use with a cleaning cloth. If necessary, use a slightly damp cloth. Ensure that no water gets inside the device while cleaning it.
- Do not use any cleaning agents or solvents, as they could damage the device.
- Under no circumstances should you hold the device under water, as this can cause liquid to enter and damage the device.

9. Accessories and replacement parts

For replacement parts and additional items visit www.beurer.com or contact the corresponding service address in your country (according to the service address list).

Replacement parts and additional items are also available from retailers.

Item	Item number and/or order number
Advanced Hydro Lift Cream (30 ml)	163.392
Intensive Cleansing Milk (200 ml)	163.394

10. What if there are problems?

Problem	Cause	Solution
There is nothing shown in the device display.	The mains adapter is not plugged in correctly.	Check that the mains adapter is correctly plugged in.
	The device is not switched on.	Press the ON/OFF button.
The heat and cool function is not working.	The mains cable is damaged.	Contact your retailer.
	The device is damaged	

11. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

12. Technical specifications

Model no.	FC 90
Input	AC 100–240 V / 50–60 Hz / 0.5 A
Output	DC 6.5 V / 2.3 A / 15.0 W
Dimensions	18.5 x 5.5 x 5.3 cm
Weight	Hand-held device: approx. 173 g Mains adapter: approx. 90 g
Average active efficiency	≥ 85.99 %
Efficiency low load	81.05 %
No-load power consumption	≤ 0.047 W

13. Warranty / Service

Beurer GmbH, Söflinger Straße 218, 89077 Ulm, Germany (hereinafter referred to as “Beurer”) provides a warranty for this product, subject to the requirements below and to the extent described as follows.

The warranty conditions below shall not affect the seller’s statutory warranty obligations which ensue from the sales agreement with the buyer.

The warranty shall apply without prejudice to any mandatory statutory provisions on liability.

Beurer guarantees the perfect functionality and completeness of this product.

The worldwide warranty period is 3 years, commencing from the purchase of the new, unused product from the seller.

The warranty only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes in the context of domestic use.
German law shall apply.

During the warranty period, should this product prove to be incomplete or defective in functionality in accordance with the following provisions, Beurer shall carry out a repair or a replacement delivery free of charge, in accordance with these warranty conditions.

If the buyer wishes to make a warranty claim, they should approach their local retailer in the first instance: see the attached “International Service” list of service addresses.

The buyer will then receive further information about the processing of the warranty claim, e.g. where they can send the product and what documentation is required.

A warranty claim shall only be considered if the buyer can provide Beurer, or an authorised Beurer partner, with

- a copy of the invoice/purchase receipt, and
- the original product.

The following are explicitly excluded from this warranty:

- deterioration due to normal use or consumption of the product;
- accessories supplied with this product which are worn out or used up through proper use (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, light sources, attachments and nebuliser accessories);
- products that are used, cleaned, stored or maintained improperly and/or contrary to the provisions of the instructions for use, as well as products that have been opened, repaired or modified by the buyer or by a service centre not authorised by Beurer;

- damage that arises during transport between manufacturer and customer, or between service centre and customer;
- products purchased as seconds or as used goods;
- consequential damage arising from a fault in this product (however, in this case, claims may exist arising from product liability or other compulsory statutory liability provisions).

Repairs or an exchange in full do not extend the warranty period under any circumstances.

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